





## Easy Implant Locating and Directing Without Surgical Stent

# Parallel Guide KIT

- Possible to set implant location and direction with Single guide, bridge guide for easy and precise drilling
- Used for wide range of surgical cases such as Single or bridge case

### Initial Drill

- The drill only for single guide
- Preventing slip in the initial drilling by adding point function
- Adjusting depth by 6~8mm (marking)

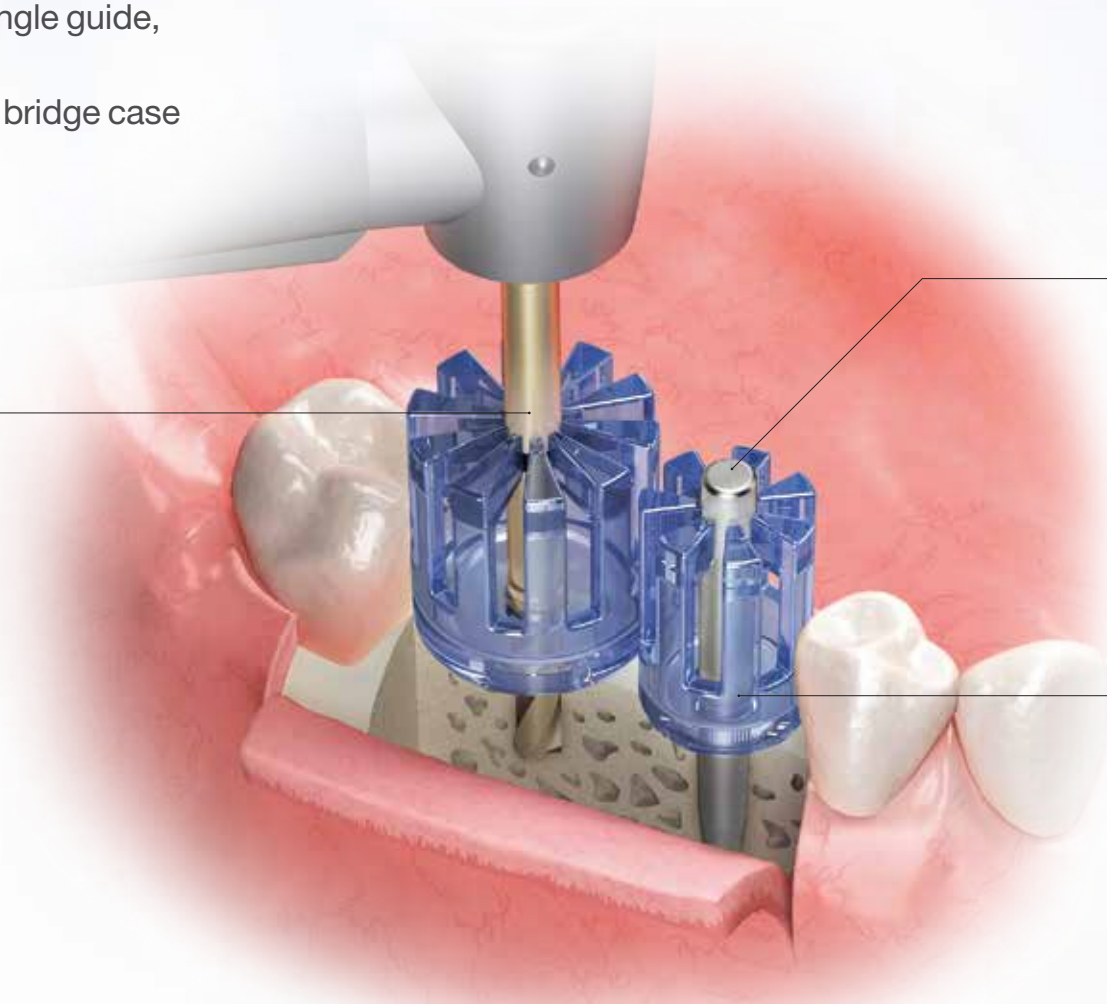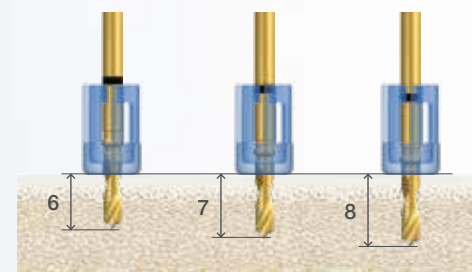

### Guide Pin

- Path checking after drilling
- Available to be used as fake teeth by connecting to single guide




### Single Guide

- Drilling guide
- Considering the size of mesio-distal crown diameters, it is designed as 1mm unit from Ø5.0 to Ø10.0
- Possible to measure the length when connected to guide pin
- No risk of infection as it is disposable

## Bridge Case

### Single Guide

- Selecting the length as the diameter of guide(diameter/2), easy to set implant location and direction by attaching to a adjacent tooth
- Available to be used as fake teeth by fastening to guide pin

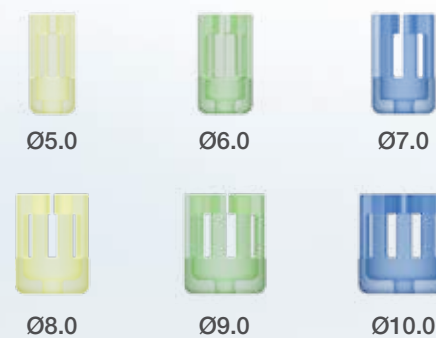

## Bridge Guide

- Using after initial drilling with single guide
- With its fan shape, easy to adjust length and direction by fixing the unit as 0.5mm from 7mm to 12.5mm

